

Urology Equipment & Accessories

For Use With HITACHI Ultrasound Systems

Micro-Touch® LP (Single-Sided Table Mount)

The Micro-Touch LP (Low Profile) is specifically designed for use on cystoscopy tables with a fixed fluoroscopy/x-ray gantry. The Micro-Touch LP attaches to tables, opposite the gantry, and allows users to perform simultaneous fluoroscopy/x-ray and seed implantation while retaining the remarkable versatility of the Micro-Touch ultrasound platform and needle guidance system.

Features:

- Rapid fixation and locking with single knob turn
- Zero gravity feel with unique weight compensator
- Precise, in vivo probe adjustment of pitch, yaw, and roll
- Intuitive design suggests a short learning curve for the clinician and staff
- Quick-connect mounting system



Micro-Touch® (Dual-Sided Table Mount)

The Micro-Touch features stable single-arm fixation with a single-point locking mechanism to instantly fix the stabilizer in the desired position. An adjustable weight-assist feature allows users to achieve a “zero gravity” feel to the entire system.

Features:

- Rapid fixation and locking with single knob turn
- Zero gravity feel with unique weight compensator
- Precise, in vivo probe adjustment of pitch, yaw, and roll
- Intuitive design suggests a short learning curve for the clinician and staff



Disposable Accessories



Sterile Drape

- Sterile, single-use disposable drape
- Latex-free, transparent material
- Cost-effective, protective barrier offers expanded sterile field



Sterile Template Grid

Disposable template grid for use with AccuCARE stepping devices.

- Sterile, single-use grid eliminates the cumbersome cleaning of reusable grids
- Easy identification and insertion of 17 or 18 GA instruments



Latex-free Endocavity Balloon

- Snug fit to transducer
- Adjustable inflation control to maximize user control over prostate position
- Assures gland positioning within the grid and at an optimal focal length to enhance ultrasound image



EXII Stepper

The EXII Stepper features a top-loading cradle and an adjustable stepper positioning reference scale. The EXII Stepper has a digital backlit monitor with an easy-to-read display which tracks the position of the transducer during the procedure. It features a serial data output port for connecting to the VeriSeed planning and treatment software.

Features:


- Incremental X and Y fine adjustment of grid to achieve precise needle path control and verification
- Allows user to set the visualized zero retraction plane at zero on the reference scale
- Selection of continuous or 5mm step positions
- Works with Varian's Variseed planning and treatment software

Classic Stepper

The Classic Stepper is a lightweight precision stepping device. The Stepper's modular design is adaptable to a variety of transrectal transducers. The Classic Stepper features a secure centerline detent featuring 45 degree clockwise and counter-clockwise rotation with easy-to-read marking scales from either side.

Features:

- Precise needle path control and verification using fine adjustment feature
- Transducer rotation out of the Stepper without loss of orientation
- Grid independently moves from a convenient backside location

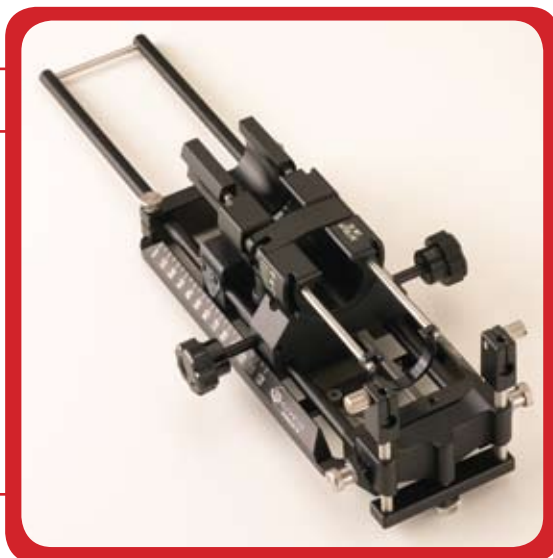

Product Ordering Information:

Steppers:

REF. #	Product Description	Qty.
644-064	Classic Stepper- Hitachi EUP-U533	1
644-065	EXII Stepper- Hitachi EUP-U533	1

Stabilizers:

REF. #	Product Description	Qty.
610-911	Micro-Touch (Dual-Sided Table Mount)	1
610-912	Micro-Touch (Single-Sided Table Mount - Right Side)	1
610-913	Micro-Touch (Single-Sided Table Mount - Left Side)	1

Disposable Accessories:

REF. #	Product Description	Qty.
610-870	Sterile Drape	20
610-898	Latex-Free Endocavity Balloon	10

Template Grids:

REF. #	Product Description	Qty.
610-905	17 GA Grid	5
610-906	18 GA Grid	5
610-909	17 GA Grid, Drape, NeoGuard Cover (2)	5
610-910	18 GA Grid, Drape, NeoGuard Cover (2)	5



Corporate Headquarters: 102 First Street South, Kalona, Iowa 52247-9589 USA
319.656.4447 | 800.445.6741 | Fax: 391.656.4451

USA Office: 1401 8th Street, Orange City, Iowa 51041 USA
712.737.8688 | 800.842.8688 | Fax: 712.737.8654

European Office: Pasteurstraat 6, 2811 DX Reeuwijk, The Netherlands
Tel. +31(0) 182-394495 | Fax: +31(0) 182-395014

COPYRIGHT © 2008. CIVCO IS A REGISTERED TRADEMARK OF CIVCO MEDICAL SOLUTIONS. PRINTED IN USA. 2008R-2328