

GOLD MARKER PLACEMENT WORKSHEET: PROSTATE

SOFT TISSUE MARKERS

The surface of our soft tissue gold marker is specially knurled, intended to inhibit migration after implantation. The markers are implanted under ultrasound or CT guidance using special needles. Implant needles come pre-loaded with one marker each, packaged and sterilized, ready to use upon receipt.



PATIENT NAME: _____
 RADIATION ONCOLOGIST: _____
 UROLOGIST: _____
 DATE: _____ NOTES: _____

Requirements:

- Sterile package containing (3) implant needles loaded with (3) one each non-migrating gold markers and plugged with sterile bone wax.
- Tweezers – clean not sterile.


Placement at base just right of urethra, midgland at extreme left side and apex slightly right as shown.

Record distance from apex and base to markers 1 and 3:

Marker 1 = _____ mm

Marker 3 = _____ mm

Lot # _____

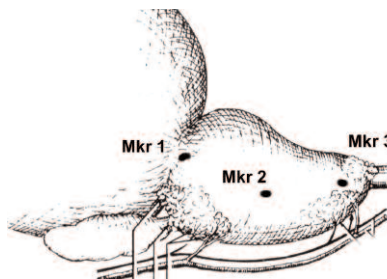

CT -MR Fusion Notes Z values

	CT	MR
#1	_____	_____
#2	_____	_____
#3	_____	_____